

磁珠純化系統


Magou

Hold住蛋白質及抗體

beads

Magnetic beads for purification and separation/isolation of bio-molecules

超方便、無需離心

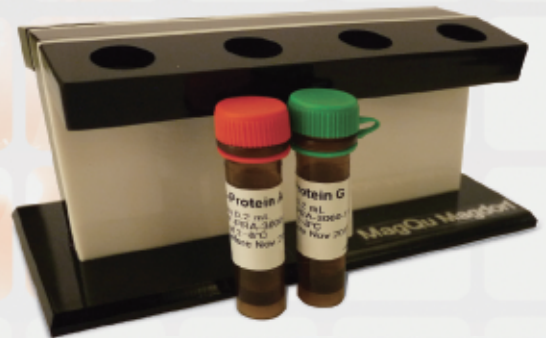

產品特色

顆粒高一致性: ~3um in diameter

核心安定性: Iron oxide Fe₃O₄


高水溶性及生物相容性: dextran-coated surface

持久可重複使用: covalent coupling of surface probe

用量極省: <5ul of beads / assay



Ts Su W1 W2 E

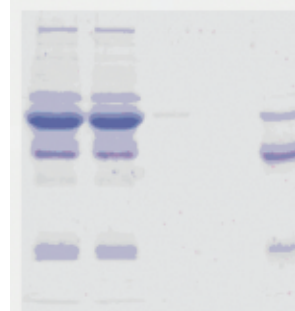

←Albumin
←Heavy chain
←Light chain

Ts : Total human serum
Su : Supernatant
W1 : 1st wash
W2 : 2nd wash
E : Elution

應用範圍

1. 抗體純化: Ig purification
2. IP and Co-IP
3. Immunoassay
4. Pull-down assay
5. Purification of His-tagged and GST-tagged protein

訂購資訊

Qbeads	Cat #	Specifications	包裝	特價
Qbeads-PRA	MF-PRA-3000	Protein A conjugate	1 mL	\$4,000
Qbeads-PRG	MF-PRG-3000	Protein G conjugate	1 mL	\$4,000
Qbeads-STA	MF-STA-3000	Streptavidin conjugate	1 mL	\$5,600
Magdorf	8管磁座	8 wells for 1.5 ml microcentrifuge tube		\$4,900
Magtractor	多孔盤磁座	Magnetic plate for 6, 24 or 96 plates		\$9,000

